



The Softest Micrus® Microcoil Ever

DELТАPLUSH™ Microcoils combine our exclusive DELTA WIND™ technology with the softest platinum wire that we have ever used. The result is a 10 system microcoil with the flexibility to find empty spaces and the softness to thoroughly fill—making DELТАPLUSH the ideal finishing microcoil.


DELTA™ Precision

The revolutionary DELTA WIND technology features a unique rounded triangular-shaped wind that provides hundreds of natural micro-deflection points, allowing the microcoil to change direction at will. The result is a microcoil with amazing softness, stability, and dense packing ability.

- Engineered with the softest Micrus .0013" wire ever used, DELТАPLUSH provides exceptional microcoil conformability
- Provides focused filling at the aneurysm neck

Finish Soft with DELТАPLUSH

- DELТАPLUSH measures 21% softer than HyperSoft*
- DELТАPLUSH measures 40% softer than Ultrasoft*

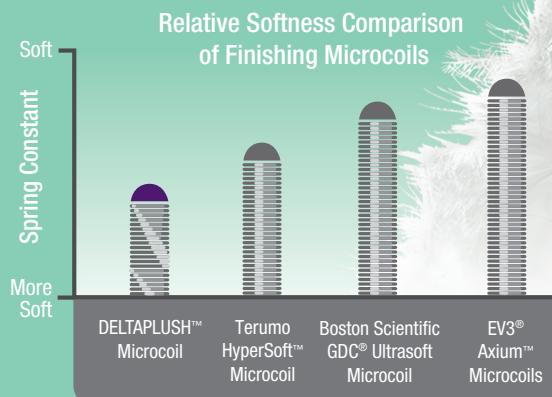


*Data on file, Micrus Endovascular Corporation.



Finish with Confidence

- 1.5 mm size provides increased treatment options
- DELTA WIND's natural micro-deflection points may reduce the tendency to coin stack while helping facilitate easier microcoil repositioning
- True 10 system outer diameter size with exceptional softness provides the ability to fill small interstices

Coil Outer Diameter (O.D.)

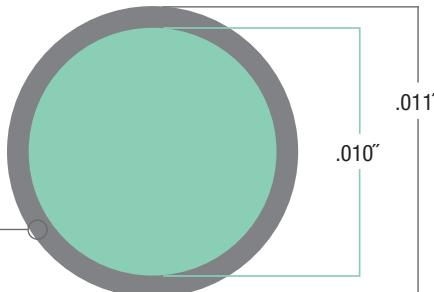

Compared with 11 system finishing coils, the coil OD of DELTAPLUSH is 10% smaller in size, helping the microcoil access the tightest spaces.

Clinical Case Experience


CASHMERE® Frame-DELTAPAQ® Fill-DELTAPLUSH™ Finish



Left ICA terminus, measuring 5.4 mm x 4 mm x 2.7 mm.



An initial 5 mm x 12 cm CASHMERE® Microcoil frames the aneurysm. A total of 18 cm of DELTAPAQ® and 10 cm of DELTAPLUSH Microcoils are used to occlude the aneurysm. Used in conjunction with Micrus framing coils, the DELTAPAQ and DELTAPLUSH Microcoils were designed to help increase overall packing density.

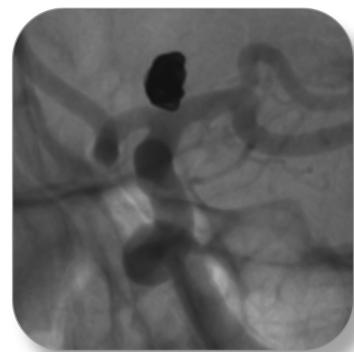

Images courtesy of Dr. S. Moskowitz and Dr. V. Gonugunta; Cleveland Clinic Foundation, Cleveland, OH



micrus[®]

endovascular

821 Fox Lane
San Jose, CA 95131
United States
T 408.433.1400
F 408.433.1401

3260 Executive Way
Miramar, FL 33025
United States
T 954.441.3000
F 954.441.3001

Avenue des Baumettes 7
1020 Renens
Switzerland
T +41.21.631.3434
F +41.21.631.3444

1st Floor, Unit 4
Stansted Courtyard
Takeley
Bishop's Stortford
Essex
CM22 6PU
United Kingdom
T +44.800.612.6363
F +44.800.612.8484

www.micruscorp.com

CAUTION: Federal (U.S.A.) law restricts this device to sale by or on the order of a physician.
DELTAPAQ, CASHMERE, MICRUS, MICRUS ENDOVASCULAR, and the M logo are registered trademarks of Micrus Endovascular Corporation.
DELTAPLUSH, DELTA WIND, DELTA are trademarks of Micrus Endovascular Corporation. Other product names or marks are the property of their respective owners.
© 2009 Micrus Endovascular Corporation. All rights reserved.

